

lionelo

User manual	Használati utasítás
Instrukcja obsługi	Manual utilizare
Bedienungsanleitung	Bruksanvisning
Инструкция по эксплуатации	Bruksanvisning
Manuale d'uso	Brugsvejledning
Manuel de l'Utilisateur	Käyttöopas
Manual de usuario	Manual de instruções
Handleiding	Οδηγίες χρήσης
Naudojimo instrukcija	Kasutusjuhend
Návod k obsluze	Lietošanas instrukcija



Milky 2in1

Milk collector 2 in 1
Kolektor mleka 2 w 1
Milchauffänger 2 in 1
Сборщик молока 2 в 1
Raccoglitore di latte 2 in 1
Collecteur de lait 2 en 1
Colector de leche 2 en 1
Melkverzamelaar 2 in 1
Pieno surinktuvas 2 viename
Sbërač mléka 2 v 1
Tejgyűjtő 2 az 1-ben
Colector de lapte 2 în 1
Mjölkuppsamlare 2 i 1
Melkesamler 2 i 1
Mælkeopsamler 2 i 1
Maitokeräjä 2 in 1
Coletor de leite com função de bomba tira leite
Συλλέκτης γάλακτος με λειτουργία θηλάστρου
Rinnarumba funktsiooniga piimakoguja
Piena savācējs ar krūts sūkņa funkciju



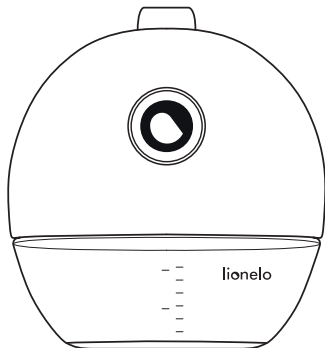
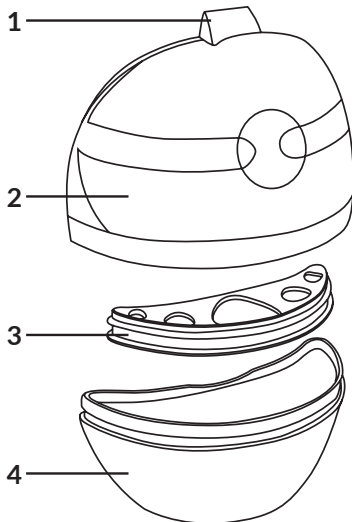
A

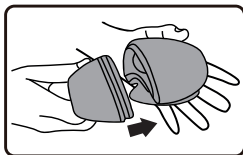
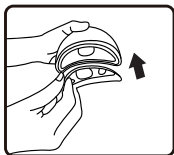
1

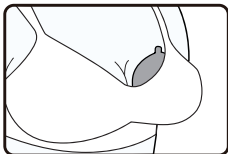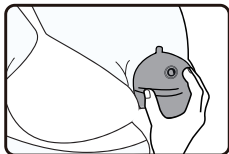2

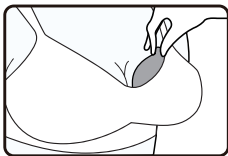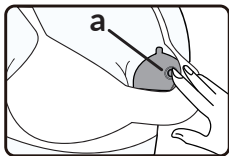3

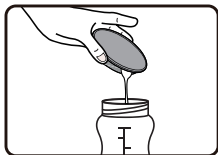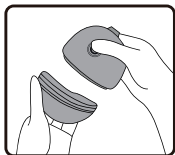4

Dear Customer!

If you have any comments or questions about the product you have purchased, please contact us:
help@lionelo.com

Manufacturer:

BrandLine Group Sp. z o.o.

A. Kreglewskiego 1, 61-248 Poznań, Poland

IMPORTANT INFORMATION

1. Before and after each use, wash and disinfect all parts that will come into contact with the body and milk.
2. For hygiene purposes, the product should be used by only one person.
3. Before each use, inspect the product for any signs of damage or defects. If any are found, do not use it.
4. While using the product, ensure the milk catcher does not press against the milk ducts. For optimal use, when positioning the funnel on your breast, ensure the nipple is centred within the funnel opening.

DESCRIPTION OF PARTS (FIG. A)

1. Airflow control opening
2. Main body with pump
3. Milk separator
4. Milk container

USE

Important! Before use, make sure the milk catcher with breast pump function is correctly assembled.

1. Insert the milk separator into the main body and then fit the milk container (Fig. 1).
2. Put on your bra, then securely fit the milk catcher with breast pump function to the breast (Fig. 2).
3. Press the button (a, Fig. 3) to express milk. The suction level of the pump can be adjusted by pressing the airflow control opening (1, Fig. A).

4. When the container is full of milk, open the airflow control opening (1, Fig. A), and then gently remove the device from the breast.
5. Detach the milk container from the main body and transfer the milk into the bottle or milk bag for storage (Fig. 4).

CLEANING AND DISINFECTION

1. Immediately after use, wash all parts that come into contact with the body and milk separately to prevent milk residues from drying out and to prevent bacterial growth.
2. To disinfect, place all parts in boiling water for at least 10 minutes. Ensure all parts are fully submerged in the water and not in direct contact with the bottom of the pot.
3. Leave all product parts to dry on a clean towel or paper towel.
4. Store the product in a clean bag, container, or other clean location until the next use. If the parts are still wet, do not store them in a sealed container.

The photos are for reference only, the actual appearance of products may differ from the photos.

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