lionelo


User manual
Instrukcja obsługi
Bedienungsanleitung
Инструкция по эксплуатации
Manuale d'uso
Manuel de l'Utilisateur
Manual de usuario
Handleiding
Naudojimo instrukcija
Návod k obsluze

Használati utasítás
Manual utilizare
Bruksanvisning
Bruksanvisning
Brugsvejledning
Käyttöopas
Manual de instruções
Οδηγίες χρήσης
Kasutusjuhend
Lietošanas instrukcija



Milky Plus

Milk collector
Kolektor mleka
Milchauffänger
Сборщик молока
Raccogliatore di latte
Collecteur de lait
Colector de leche
Melkverzamelaar
Pieno surinktuvas
Sběrač mléka
Tejgyűjtő
Colector de lapte
Mjölkuppsamlare
Melkesamler
Mælkeopsamler
Maitokerääjä


Coletor de leite
Συλλέκτης γάλακτος
Piimakoguja
Piena savācējs

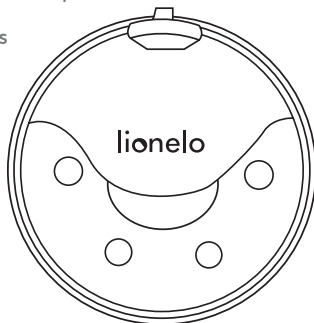
www.lionelo.com

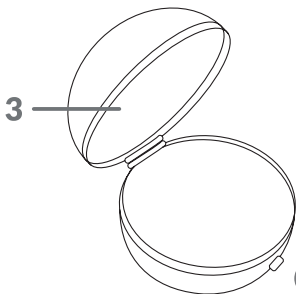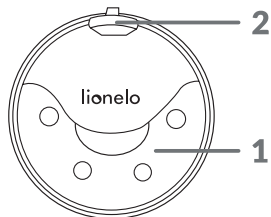A

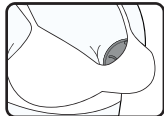1

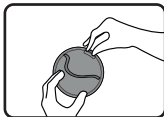2

Dear Customer!

If you have any comments or questions about the product you have purchased, please contact us:

help@lionelo.com

Manufacturer:

BrandLine Group Sp. z o.o.

A. Kreglewskiego 1, 61-248 Poznań, Poland

IMPORTANT INFORMATION

1. Before and after each use, wash and disinfect all parts that will come into contact with the body and milk.
2. For hygiene purposes, the product should be used by only one person.
3. Before each use, inspect the product for any signs of damage or defects. If any are found, do not use it.
4. While using the product, ensure the milk catcher does not press against the milk ducts. For optimal use, when positioning the milk catcher on your breast, ensure the nipple is centred within the opening.

DESCRIPTION OF PARTS (FIG. A)

1. Catcher x2
2. Silicone pin
3. Storage case

USE

1. Put on your bra, then securely fit the milk catcher to the breast (Fig. 1).
2. Once the catcher is full of milk, carefully remove it from the bra. Detach the silicone pin and transfer the milk into a bottle or milk bag and place it in the fridge (Fig. 2).

3. Afterwards, wash the catcher thoroughly and allow it to air dry.

CLEANING AND DISINFECTION

1. Immediately after use, wash all parts that come into contact with the body and milk separately to prevent milk residues from drying out and to prevent bacterial growth.
2. To disinfect, place all parts in boiling water for at least 10 minutes. Ensure all parts are fully submerged in the water and not in direct contact with the bottom of the pot
3. After washing, leave all product parts to dry on a clean towel or paper towel.
4. Store the product in the attached container or other clean location until the next use. If the parts are still wet, do not store them in a sealed container.

The photos are for reference only, the actual appearance of products may differ from the photos.

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